




Automated Blood Pressure Monitor
for Companion Animals



SunTech® *Vet20*TM
BP MONITOR

A Reliable Vet BP Monitor That is Fast, Simple & Silent





Two companion animal modes

Memory   

#	Sys	Dia	HR	MAP	Animal
10	122	80	110	94	Large
9	118	79	85	95	Large
8	 Air leak in hose				
7	119	80	95	94	Large
6	120	77	100	95	Small

Up to 50 stored data sets

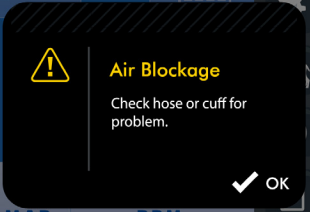

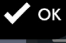Air Blockage
Check hose or cuff for problem.


MAP BPM

Troubleshooting guide



+ Motion-tolerant technology
equals fewer failed readings


+ Animal-specific blood pressure
algorithm

The new portable **SunTech® Vet20™** quickly and easily provides reliable BP measurements on awake companion animals, ensuring a simpler assessment process for veterinarians and their staff members. Utilizing *SunTech's* trusted Advantage™ VET BP technology, the *SunTech Vet20* quietly performs motion-tolerant, oscillometric BP measurements, delivering more accurate results with a greater success rate than other devices. Users are able to get started using the device immediately as the intuitive touch screen interface includes all of the necessary features for reviewing and analyzing BP data, yet little to no staff training is required. Alarms can be silenced to avoid frightening the animal and saved measurement data is always retained when the device is powered off.

+ Features & Benefits

+ Description

- Reliable:** As a result of *SunTech's* motion-tolerant technology, the *SunTech Vet20* produces more successful BP readings than other devices.
- Simple to Use:** Intuitive, touch screen interface has all necessary features for reviewing and analyzing BP data and requires little to no staff training.
- Accurate:** BP Algorithm was specifically developed for companion animals by experts at a leading vet school.
- Silent:** Device quietly measures BP and alarms may be silenced to avoid frightening the animal.
- Fast:** Measurements are done faster than competition. Doppler takes 10 minutes!
- Portable:** Battery-powered device can be easily moved between exam rooms.

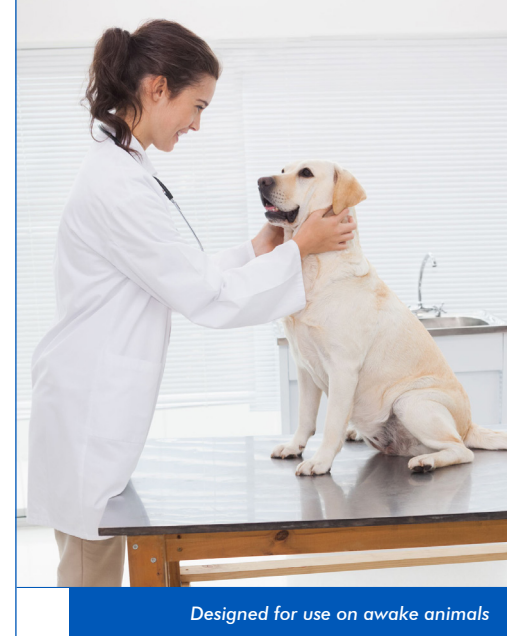

and there's more...

The *SunTech Vet20's* **advanced memory** ensures data is always retained when the device is powered off. Users can also select specific measurements from the device memory and perform **automated BP averaging**.

+ Specifications

+ Description

+ Cuffs

- | | | |
|-----------------------|----------------------------|--------------------|
| SYS: | 40 to 265 mmHg | #1 3-6 cm |
| DIA: | 20 to 200 mmHg | #2 4-8 cm |
| MAP: | 27 to 222 mmHg | #3 6-11 cm |
| HR: | 25 to 300 BPM | #4 7-13 cm |
| Device Weight: | 1.3 lbs (with batteries) | #5 8-15 cm |
| | 590 grams (with batteries) | #6 12-19 cm |
| Dimensions: | 6.25" x 5" x 5.25" | #7 17-25 cm |
| | 15.9 cm x 12.7cm x 13.3 cm | |
| Warranty: | 2 years on monitor | |

Softer, latex-free materials and rounded corners provide a gentle cuff solution. SunTech Medical BP cuffs are color-coded for quick size selection and range indicators ensure accurate cuff sizing and fit.



SunTech Medical, Inc.
507 Airport Boulevard, Suite 117
Morrisville, NC 27560-8200 USA
Tel: + 1.919.654.2300
1.800.421.8626
Fax: + 1.919.654.2301

SunTech Medical, Ltd.
Oakfield Industrial Estate
Stanton Harcourt Road
Eynsham, Oxfordshire OX29 4TS
England
Tel: + 44 (0) 1865.884.234
Fax: + 44 (0) 1865.884.235



www.SunTechMed.com

SunTech Medical, (Shenzhen) Co., Ltd.
2/F of Building A, Jinxiongda Technology Park
Guanlan, Bao'an District Shenzhen, 518110, PRC
Tel: + 86-755-2958 8810
+ 86-755-29588986 (Sales)
+ 86-755-29588665 (Service)
Fax: + 86-755-2958 8829

www.CN.SunTechMed.com

